

EM-1201

Twelve-channel ECG



EM-1201

Twelve-channel ECG



Wireless Network Adapter*



Barcode Scanner*



LAN Port*



USB Printer*



U-disk*



Keyboard*




SD Card*



Large Foldable LCD Display, Clear and Intuitive

10 inch TFT color display, which can be folded from 0-90°

Background with grid, 12 leads ECG can be displayed on the same screen

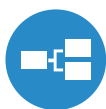

Peripherals and Communications

Integrated SD card slot and USB port for connection to barcode scanner, wireless network adapter, U-Disk, SD card, keyboard



Flexible Printing Solutions


Built-in high-resolution thermal printer, supports roll and Z-fold paper

Direct Connectivity to External Printer*


Supports ECG report network printing*

• Software Features

- Simultaneous 12 leads ECG acquisition and display
- Mapping data directly via USB cable to PC
- Supports for transmission of ECG data to ECG workstation* or ECG manager*
- long-term real-time ECG waveform freezing and review to assist diagnosis
- With paper-save function, ideal for routine ECG examination.

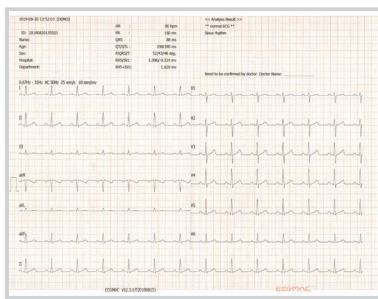


6CHx2+1R

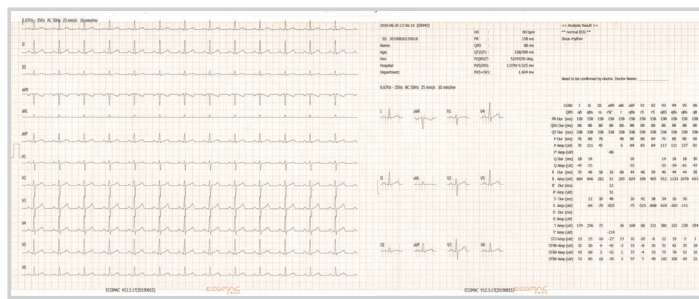


12CHx1

• Seamless Connectivity

- Multiple format reports export, such as PDF, JPG, BMP, EM-XML, SCP*, DICOM, FDA-XML*

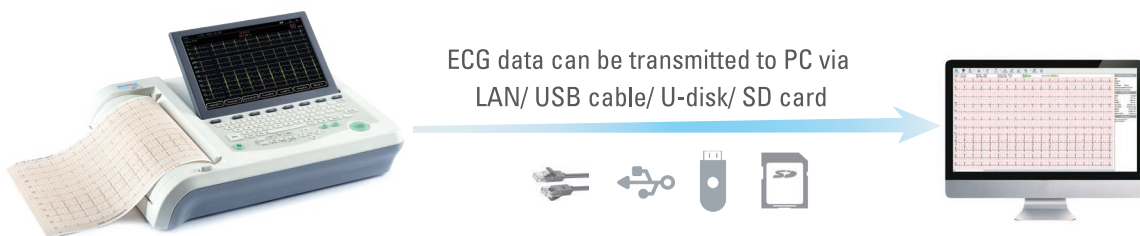

- Data access with HIS/ PACS/ EMR system based on FTP* / DICOM* / HTTP* standard protocol
- Patient information can be retrieved from the worklist*

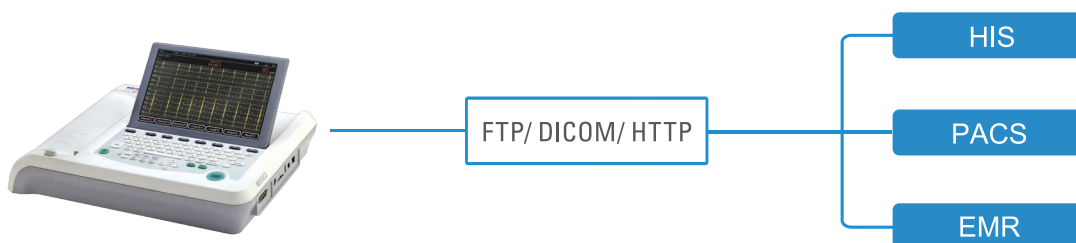

Items with * are optional

Better Care for Human Heart Health with Science and Technology
www.ecgmac.com

About us

Since 2009, Shenzhen ECGMAC Medical Electronics Co., Ltd (ECGMAC) has focused on innovating and developing cardiac diagnostic and monitoring products. As of now, we have expanded our product lines to include Multi-channel ECG, PC-ECG, Bluetooth ECG Holter and Ambulatory Blood Pressure Monitor. In the past decade, ECGMAC has 40 patents and software copyrights. ECGMAC is authenticated by German TUV ISO13485 Quality Management System (including R&D system) and Japan Quality Management System, and products have both CFDA and CE approval. With our high-quality and innovative medical devices, ECGMAC recorded a steadily increasing demand in Australia, Japan, Europe and many other countries in the world.

Shenzhen ECGMAC Medical Electronics Co.,Ltd.

Office Add: Rm 1101,1110 Baoli Building, No. 2702 Nanhai Road, Nanshan District, Shenzhen, China

Factory Add: F15-01, Block C, Konka Guangming Technology Center, No. 288 Xingxin Road,
Guangming District, Shenzhen, China

Tel: +86-755-27697821 27697823 27948579

Fax: +86-755-27697823-616

Email: info@ecgmac.com

©Shenzhen ECGMAC Medical Electronics Co.,Ltd,All Right Reserved we reserve the right to change,update,revise the product design, specification and fuction at any time without advance notice. we take no responsibility for any consequence caused by the publication